

Comfort Digisystem product information

Receiver DE10



Receiver with crystal clear sound

Receiver DE10

Small, flexible receiver which is attached directly to the ear; usable with or without all-in-the-ear-hearing aids. The receiver can be used on both right and left ear.

Suits users with high demands for nicely designed and discreet products, looking like one of those modern mobile phone accessories. The volume can be easily adjusted with a button.



Unique characteristics

SST (SecureStream Technology) – Digital sound transfer
This unique digital sound transfer between the transmitter and the receiver eliminates all transmission noise.

Moreover, the signals are encrypted, so no one can overhear what is being said.

Unlike other conventional encrypted systems, SST is free from time lag and crystal clear.



Comfort Audio
HEAR THE FUTURE

Technical facts

Receiver DE10

General

Size: 59,6 x 50,7 x 19,6 mm incl. earloop
Size: 56,7 x 23,4 x 19,6 mm excl. earloop
Weight: 15 g

Power supply

Battery: Chargeable Li-Polymer
Operating time: Up to 6 hours
Charging: Overcharging protection, automatic float charging
Charging time: Approx. 2 hours, 90% after 1 hour.

User interface

Settings: Via button
Function indicator: Easy multi-coloured status information
Manual: USB or www.comfortaudio.com

Special functions

Volume control: at the long side of the unit
Max volume: Option to set max volume; volume control is then limited

Radio Transmission

SST: Secure Stream Technology – FSK, encrypted digital sound transfer.
Frequencies: DG 03001 (F1) 818 – 872 MHz (complete package DG20001 F1)
Range: Approx. 30 m


Technical specifications

Type of receiver: SST
Antenna: Built-in



Functions



Applications:

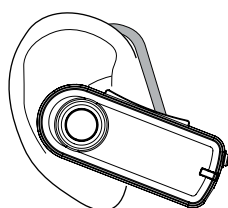

Ideal setting at the ear

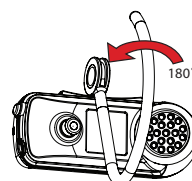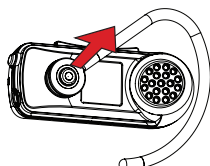

Turn around earloop

